

After-sale Services

- 1 Product free warranty period is **One year**, the original power adapter is an accessory, free warranty is **Six months**. Goggles and Mouth retractor are consumable items, free warranty is **Three months** ,Plastic LED head cover is a “one-time-use” product with no warranty.
- 2 Our company provides the warranty according to the serial number and accessory production date/lot. If the product does not have bar code label or any warranty certificate, it will be out of our free warranty;
Do NOT erase the numbers(Label)on the machine.If erased or modified, our company has the right to refuse warranty.
- 3 The following is not covered by warranty, implemented paid maintenance,if;
a, the product does not have S/N, or S/N is broken/modified
b, man-made damage, improper operation or irresistible natural disasters caused failure
c, user opens/changes/repairs the products by themselves without notifying our company or authorized dealers.
d, over the warranty period

Product Warranty Certification

Model NO.:		Product Date:	
User name:		Address :	
Tel:		E-mail :	

Skylun Medical Instrument Co.,Ltd for PrevDent International BV
Address: No.51, 1st Street Jiuliumu, ZhongCun Town, Panyu Dist.,
Guangzhou City, Guangdong Province , P.R.China.
Postal Code: 511487
Tel: 86-(20)-39204875
Fax: 86-(20)-39204875
info@skylun.com
www.skylun.com

Ver.B/3.0



Standard Accelerator

SL-TW900 User Manual



In order to reach a better result,please read this instructions before use.

Sound Warning

- a) . After starting, the machine will beep each minute.
- b) .You will hear three beeps after it has finished.

Cooling Setting

- a) . To guarantee that all LEDs work at normal temperature during whitening process, we fixed a high frequency and quiet fan to expel heat.
- b) . Built-in smart temperature protection. When high temperature is reached,it will shut down power automatically, after cooling, everything works properly.

Liability

Assembly, stretching, readjustment, modification and repair are performed by a company authorized person.
Please follow electrical products installation relevant requirements.
Please follow the instructions to use this equipment. .

Catalogue

1.Safety	2
2.Application	2
3.Installation Sketch Maps	3
4.Product Structure	6
5.Product Feature.....	8
6.Product Specification.....	8
7.Operation Instruction.....	9
8.Liability	16

Safety

Note:


Prior to installation and unit, please read the instruction carefully.
Please note the following safety use recommended in the instructions.

- ⚠ Caution: It is important to observe the following precautions and procedures.
- ⚠ Caution: Please do not use this device expose on people, who have a history of photosensitivity or photosensitizing drugs.
- ⚠ Caution: This equipment is not suitable for use in the presence of flammable anesthetic Mixtures with air or nitrous oxide.
- ⚠ Caution: Please do not look at the light directly, remember to wear goggle when it is in use, or it will hurt your eyes.
- ⚠ Caution: Please release the locking nut and handwheel when you adjust the position of the machine, and then readjust, or it will influence the life of the machine.


Application

PrevDent Standard Accelerator is a lamp delivering a specific temperature needed to accelerate the whitening agent in the PrevDent CrWR teeth-whitening procedure.

8.Screen Brightness Setting Instructions

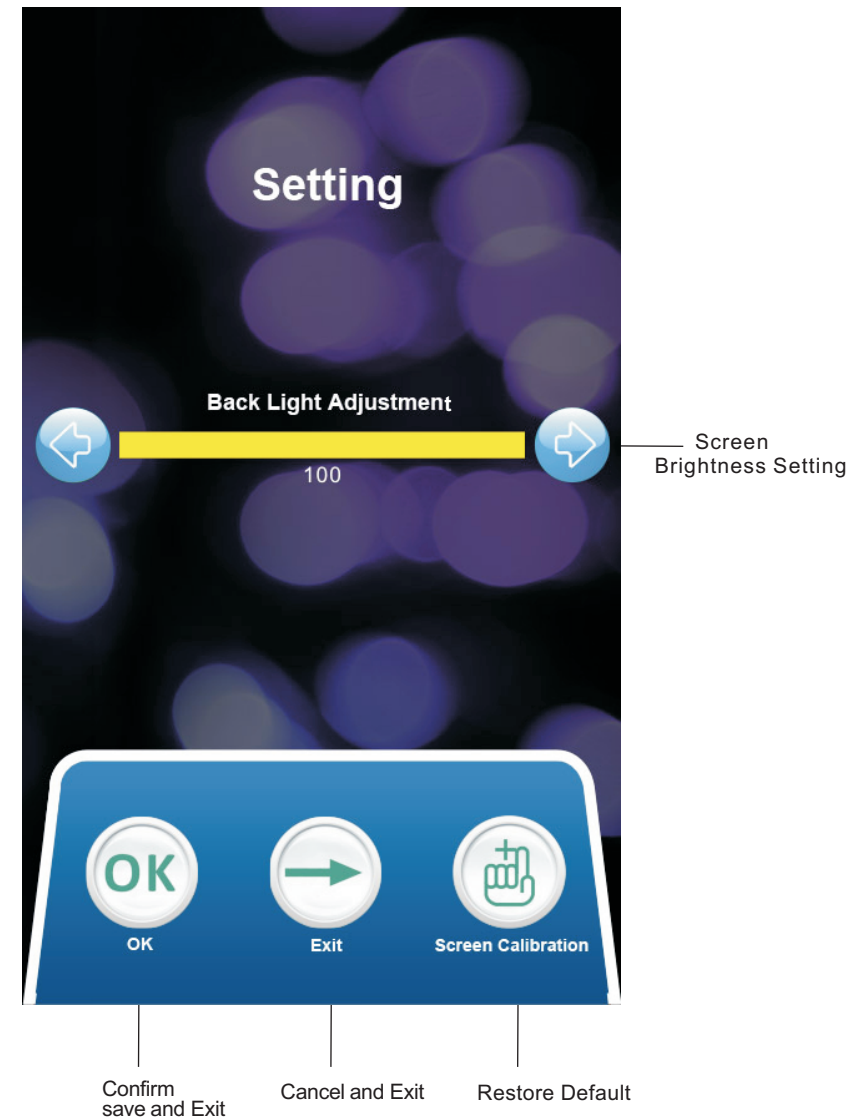

7. Mode Setting Instruction

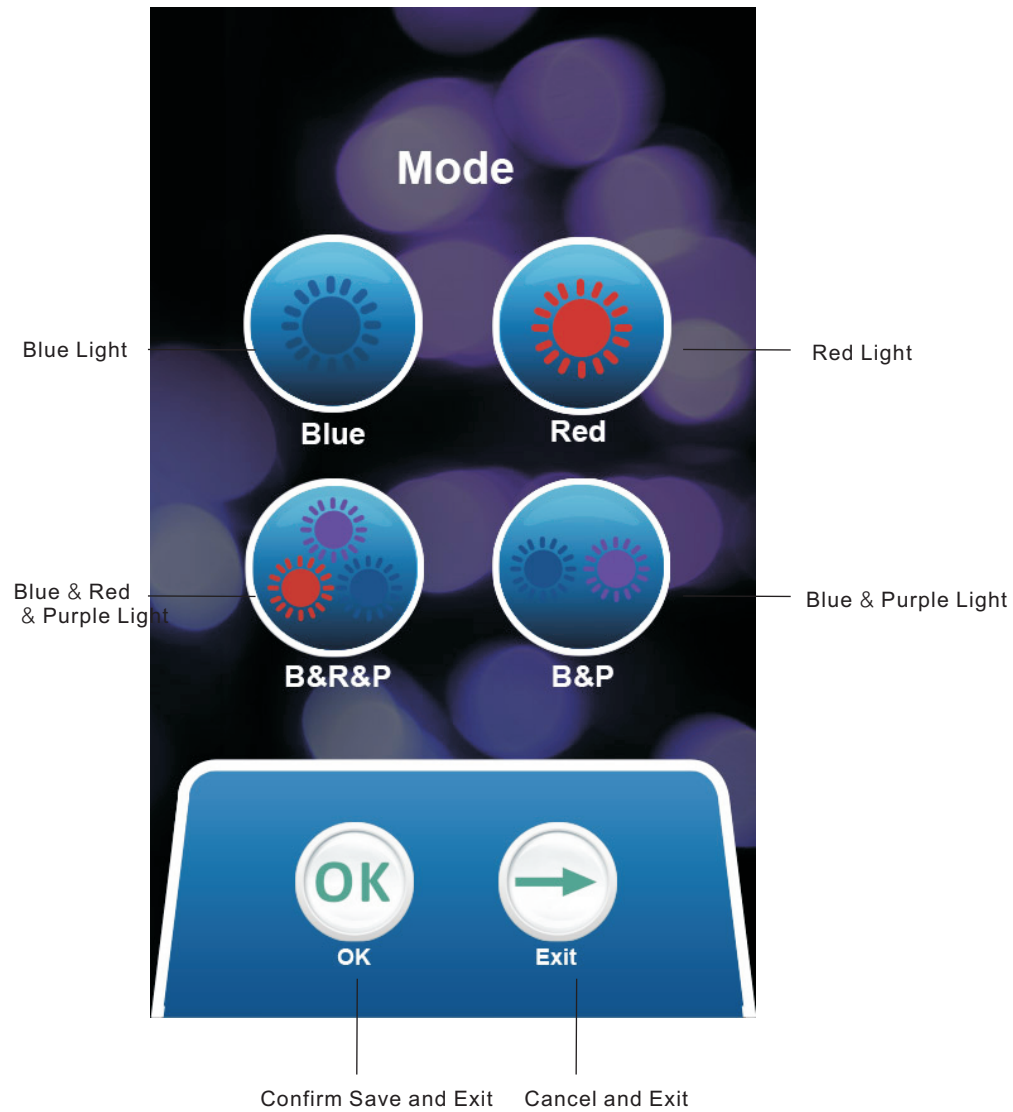

Installation Sketch Map



1. Twist the wheels to the feet, and twist tightly with a wrench

(1)



2 .Fix 5 feet onto the settled seat.

(2)


Installation Sketch Map

3. Put pole 1 with decorative set placed on the settled seat and support with hand.



(3)

6. Time Setting Instruction

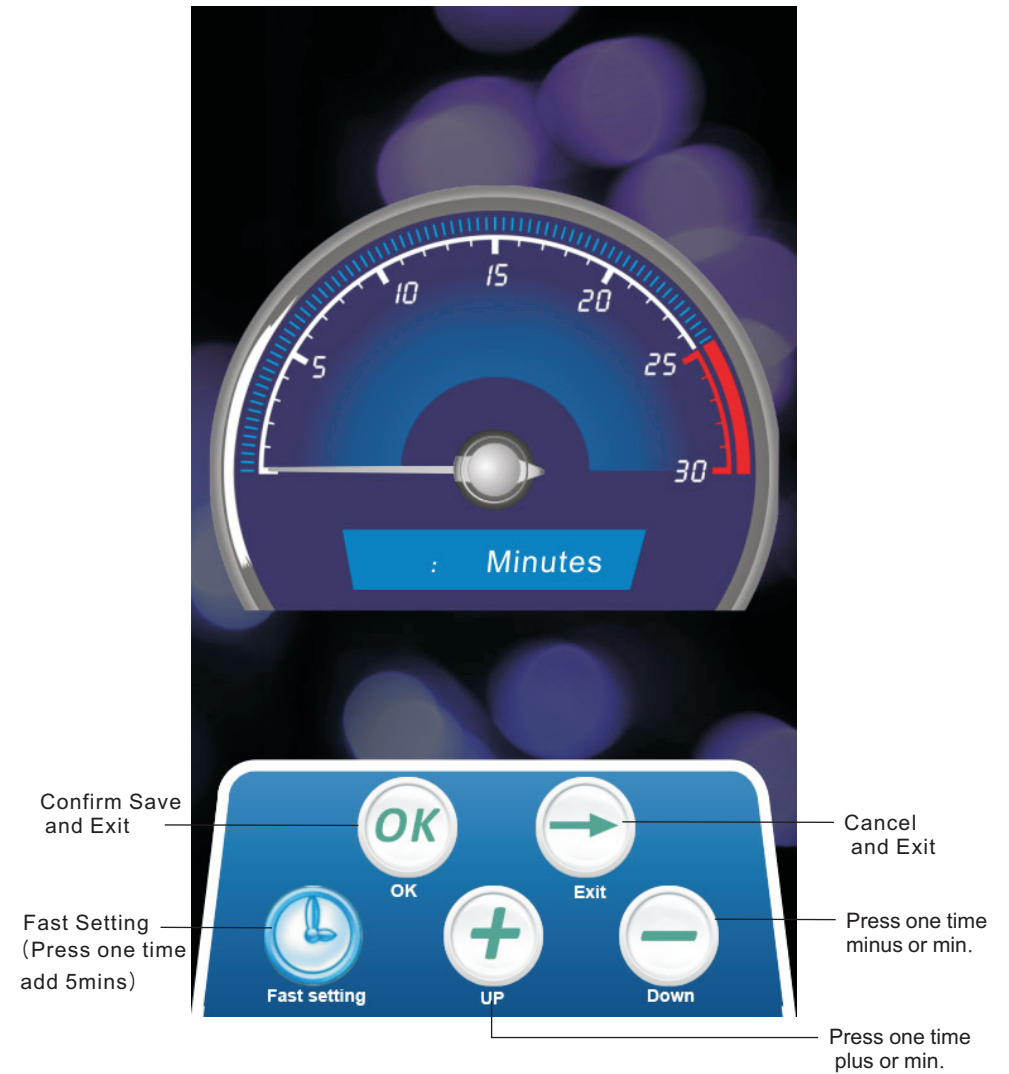

5.Power Setting Instructions (Set brightness for light output)



Installation Sketch Map



4. Turn over the feet and fix settled seat and pole with long screw.
(if the screw hasn't been tightened enough, the pole will shake)

(4)




5. Put pole 2 into pole 1 and twist nut tightly.

(5)

Product Structure

⚠ Please release the nut first, if you need to adjust machine position or angle otherwise it might effect machine using life.

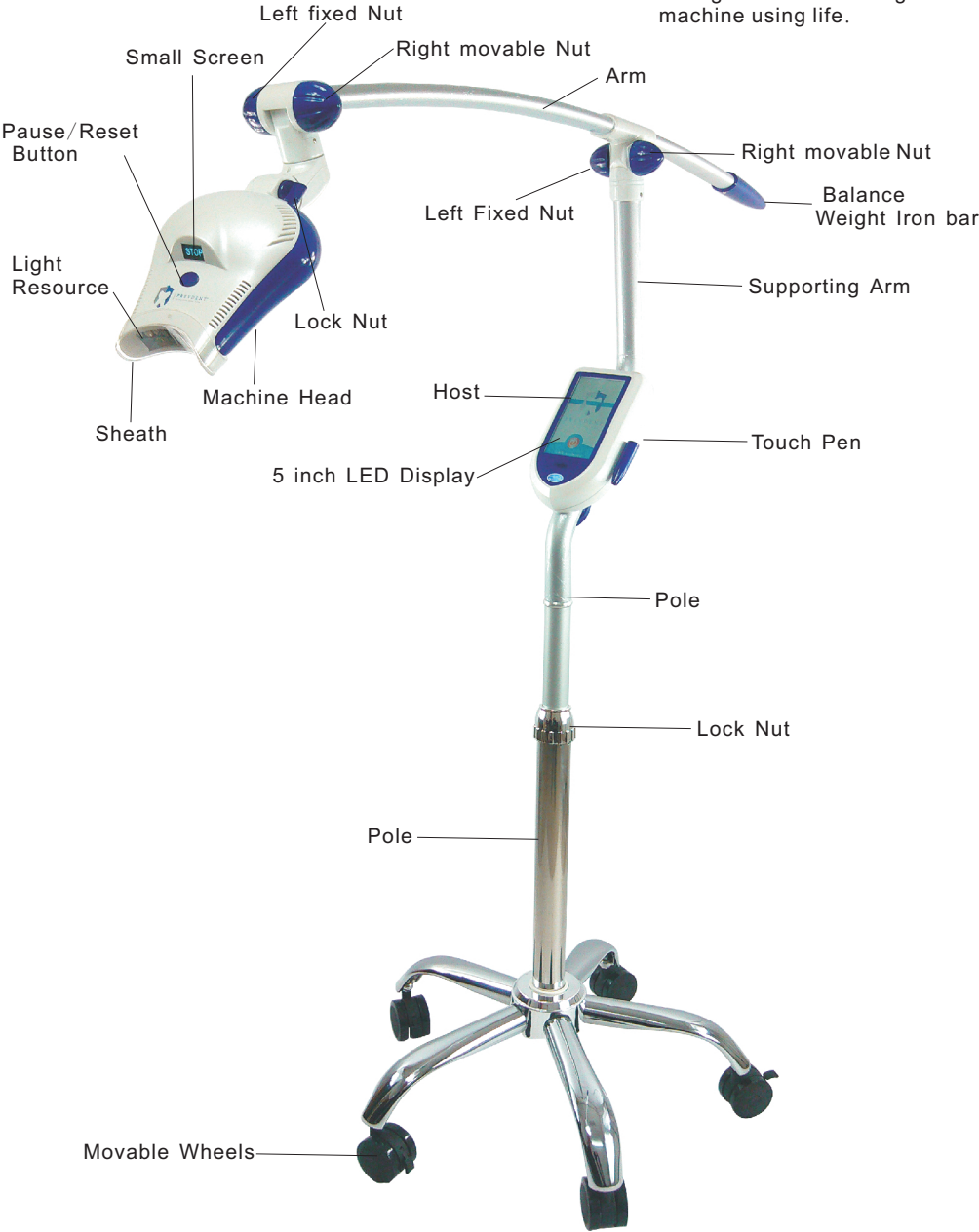


4. Main Screen Function Instruction



3. Power On Screen instruction



Host Back instruction

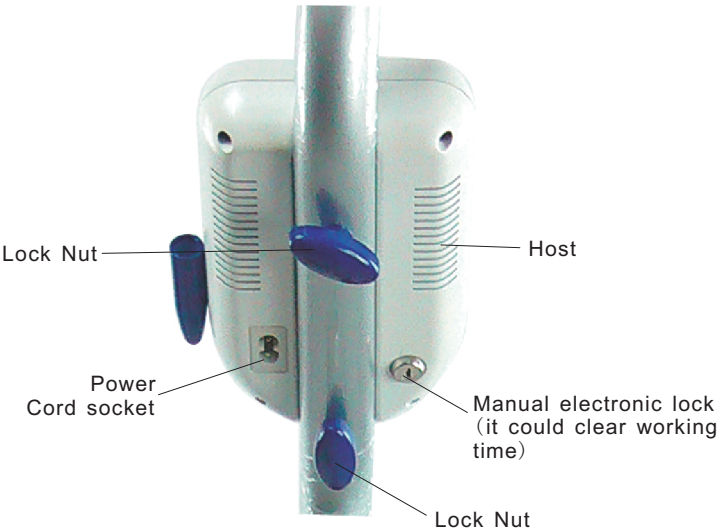

Accessories

Main unit	1pc
Power Cord	1pc(2m)
Safety Goggles	2 pairs
Sheath	20pcs
Mouth Retractor	1pc
User Manual	1pc


Product Feature

- * High efficiency wall conjoined sliding structural design, adjust angle freely and convenient to use.
- * A multi-functional high-sensitivity 5-inch touch electronic display, easy operation
- * Four light source output selection: blue, red, blue + red + purple, blue + purple
- * Three-stage power output selection, users can select a different power according to the patient's condition
- * The user can select needed time according to desired treatment, time adjustment range from 1 to 30 minutes
- * Work time counting function, comes with electronic lock, work time can be cleared by user with the lock.
- * If not used, the machine will get to sleep mode automatically after 60 seconds.
- * Enclosed disposable sheath, prevent cross-infection, clean and convenient, let patient feel relieved.

Product Specification

- * Broad Spectrum : 450nm~510nm (Blue Light) & 620nm~640nm (Red Light)
380nm~420nm (Purple Light)
- * Voltage : AC: 100V~240V
- * Light source : 4~5W/pc (Blue Light) ; 2~3W/pc (Red Light); 1W/pc (Purple Light).
- * Gross Weight : 5kgs.
- * Total Arm Length : 90cm.
- * Total Height : 120cm.

Distance between Teeth and Illuminator (mm)	Uniformity
10mm	Over 95%
20mm	Over 90%
30mm	Over 85%
40mm	Over 80%
50mm	Over 75%

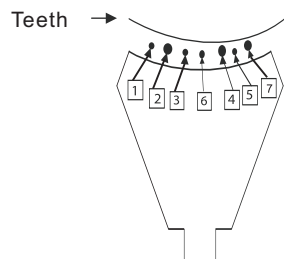


Operation Instruction

1. Turn on Teeth Whitening Machine power

- Plug power cord in.
- Turn On whitening machine power switch.



2. Head Introduction

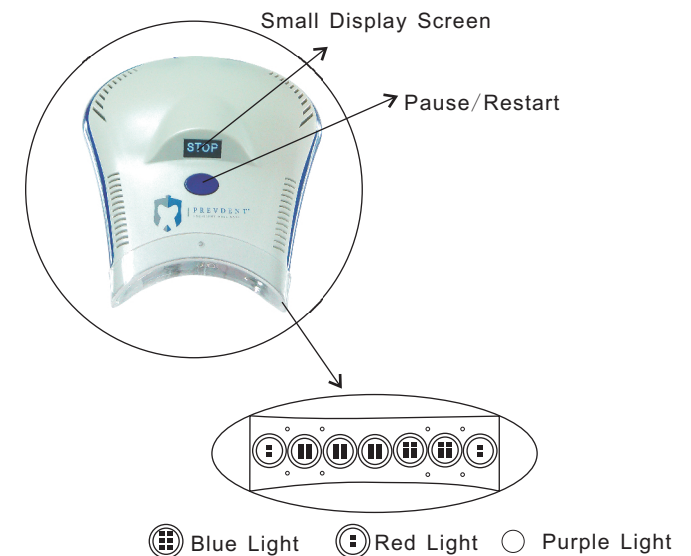

Contain five high-power Blue LEDs, two Red LEDs and 8 Purple LEDs